

Elite

Apogee 3500

Digital Color Doppler Ultrasound Imaging System

*Universal Genius
Supports
All Your Needs*



As a color Doppler ultrasound system, Apogee 3500 Elite inherits SIUI's proven imaging quality, and commits itself to be a universal genius to support all your needs. With the fashion and compact design, it features with versatile application solutions and user-friendly workflow to ensure professional diagnosis, bringing you more confidence and convenience.

Comprehensive


Equipped with various probes which self-developed by SIUI, Apogee 3500 Elite guarantees you reliable diagnosis to meet different requirements with advanced image processing technology and integrated application packages.

Advanced 4D experience in OB/GYN

The system supports both abdominal and transvaginal volumetric probes to realize 3D/4D imaging in OB/GYN application, providing more detailed and complete volumetric information to assist the observation and diagnosis of fetus wellness and uterus health.




Fetal face 4D

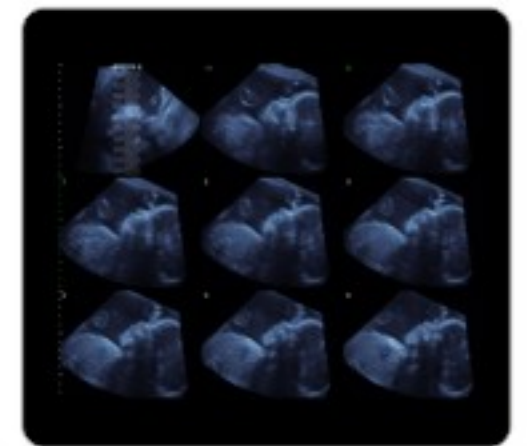


nSlice

Smart Elastography

The system supports linear probes with elastography images to visualize the stiffness of tissues in real time by delivering an external compression on the tissues. With the smart Elastography developed by SIUI, the doctors can feel more convinced in early detection of latent disease.

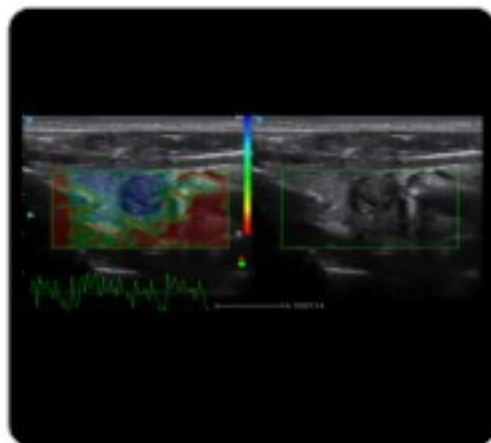


Thyroid Elastography

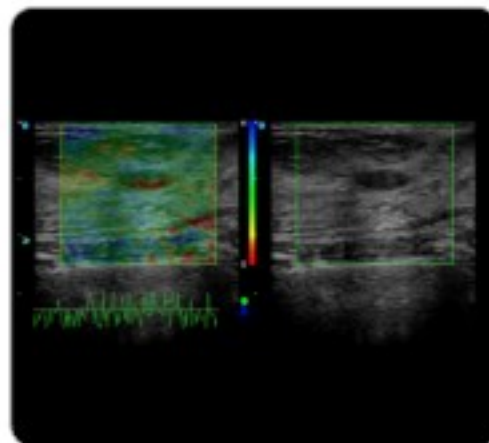


Mammary Elastography

Advanced Cardiac Functions

Including color tissue Doppler and pulse wave tissue Doppler, TDI helps to assess the directional and temporal phase of cardiac, so as to display the speed of the vessel wall and cardiac movement. Moreover, with free 360 degrees rotation and up to 3 sample lines option, the powerful anatomical M mode assists more exact analysis of cardiac structure movement even in difficult heart positioning.

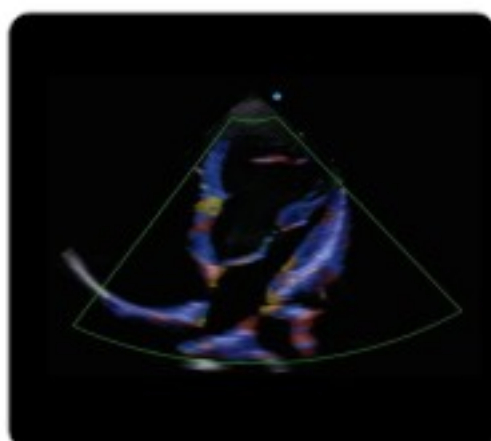

TDI

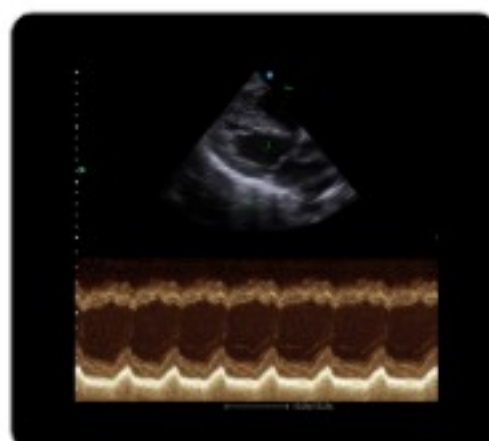

AMM




Intuitive

Touch screen operation and intelligent workflow offer concentrated and convenient diagnosis.

Auto fit

One button optimization function adjusts TGC and B gain in B mode as well as base line, PRF and PW gain to realize the best B/W image and PW image.

B mode

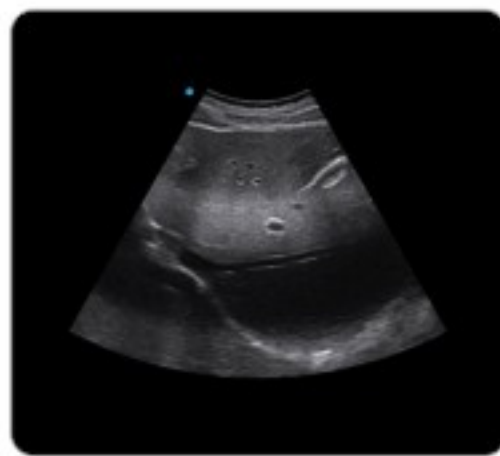


Before




After

PW mode

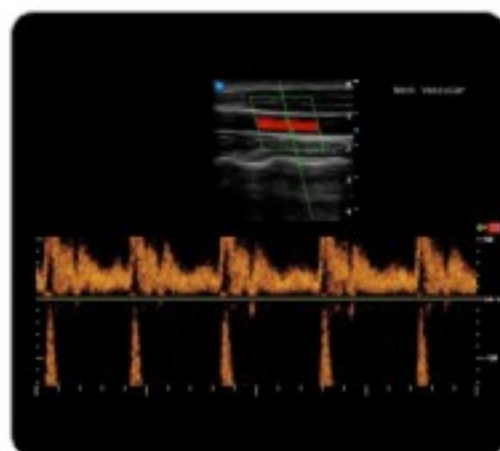



Before

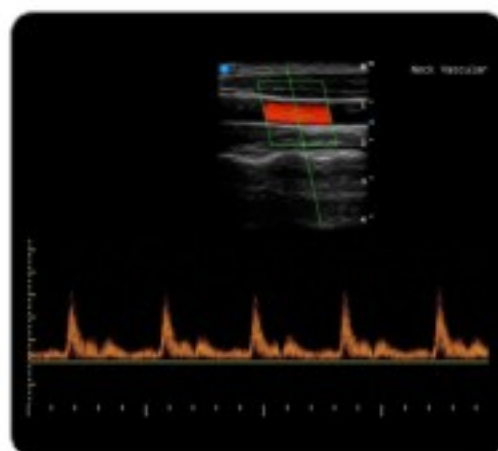


After

Auto EF

By setting 3 points of the left-ventricle (systole-diastole), system can automatically track the edge of myocardium, calculate the LV volume and output function.

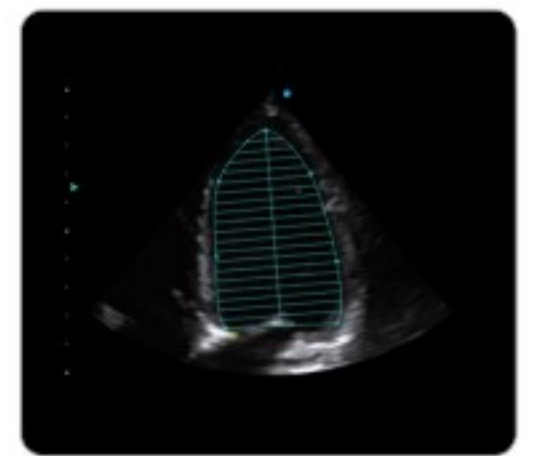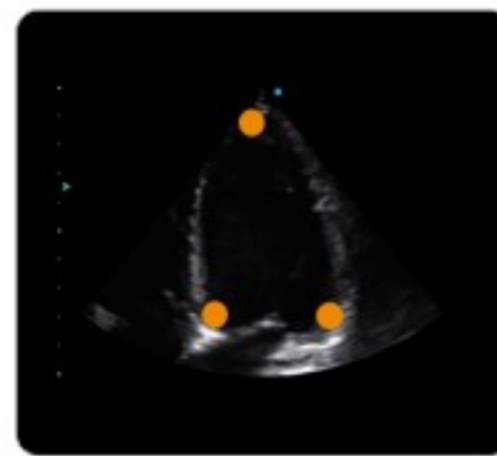

4Ch Systole

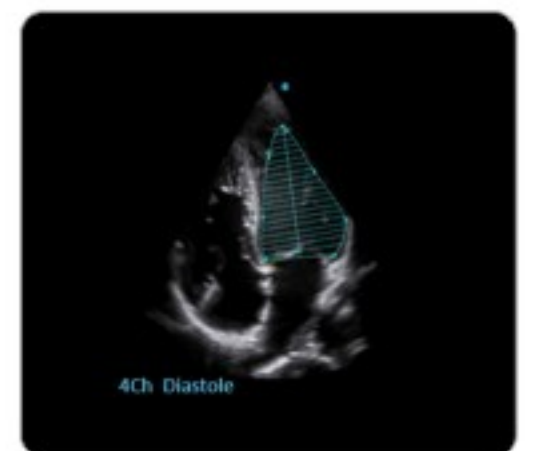

4Ch Diastole

Left-ventricle Auto EF

Auto IMT measurement

Auto IMT helps measure the Intima-Media Thickness of the carotid artery wall automatically, so as to evaluate cardiovascular disease such as hypertension diabetes.



Carotid artery Auto IMT

Smarchive

The smart patient info management system enables edit, manage and transfer patients' information as well as images, so as to realize archive, analysis and review patients' information and make remote diagnosis available.



Smarchive screenshot


Ergonomic

The design of Apogee 3500 Elite follows the idea of compact and simple to offer comfortable and pleasing operation for users.

18.5" rotatable LCD monitor



8.4" touch screen

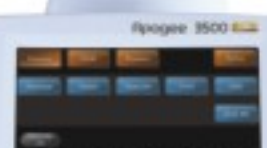

USB port



Independent probe and cable management



Swiveling control panel



Four active transducer connectors



Image Gallery


Apogee 3500 Elite is empowered with advanced image processing technology like MFI, VS-Flow, Nanoview etc. to present quality images.




Thyroid



GB Crystallization

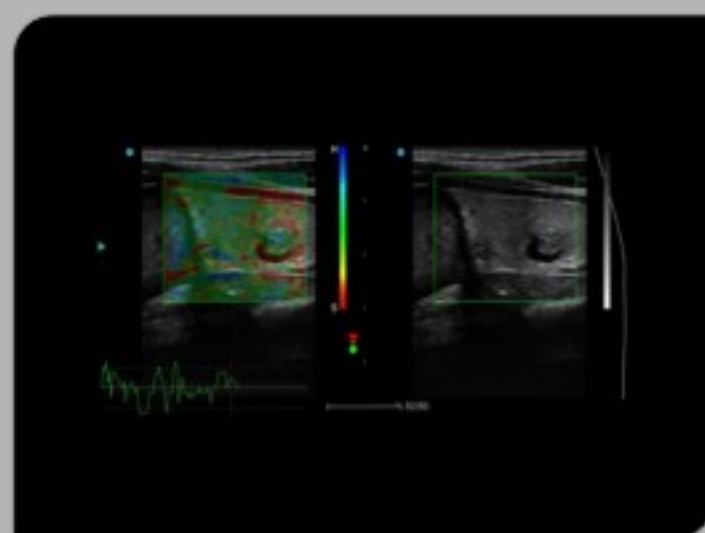


Thyroid Nodule Elastography

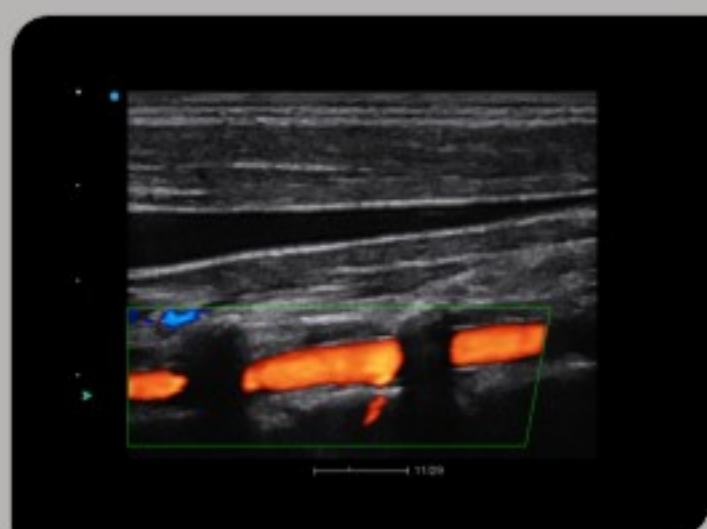


Vertebral Artery




Kidney VS-Flow

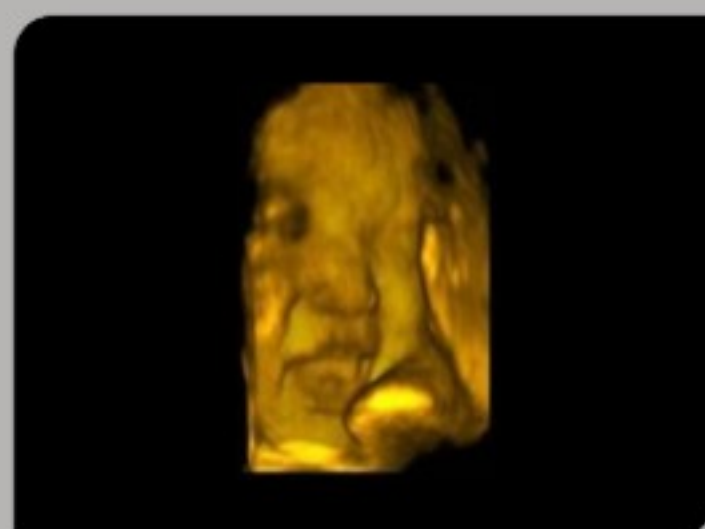


4D Body face

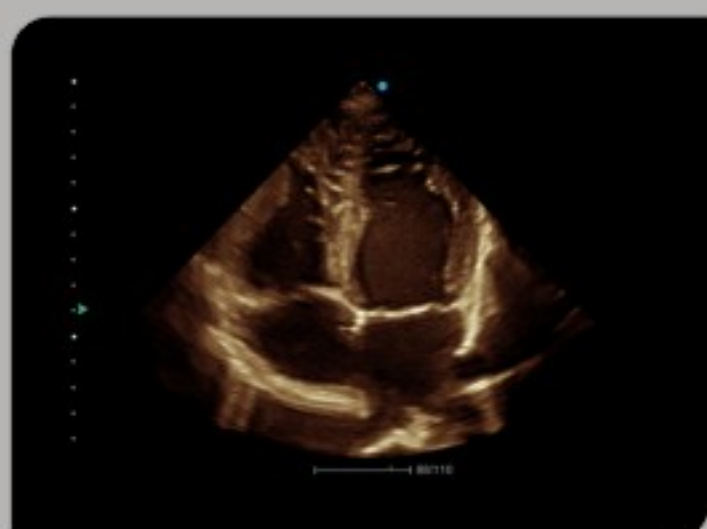

Four Chamber

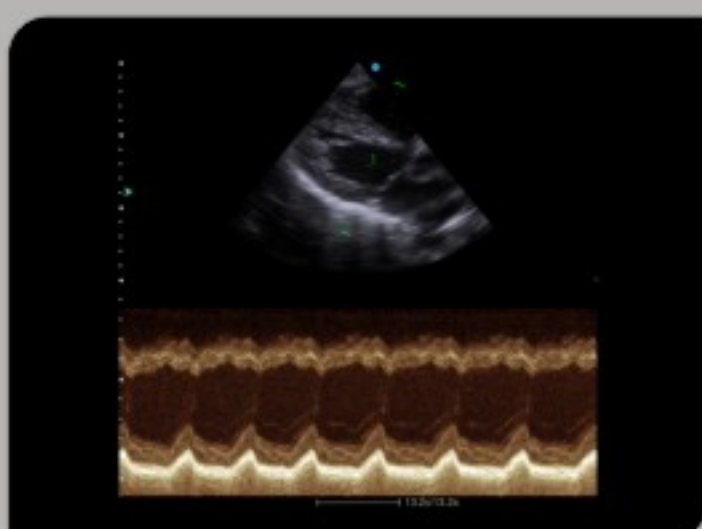

AMM

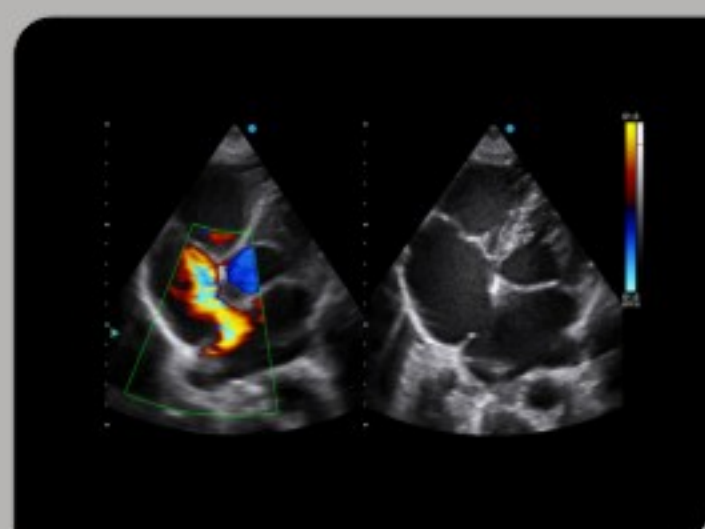

ASD



Shantou Institute of Ultrasonic Instruments Co., Ltd.

HEADQUARTERS:

Shantou Institute of Ultrasonic Instruments Co., Ltd.

Add: No.77, Jinsha Road, Shantou 515041 Guangdong, China

Tel: 86-754-8825 0150 Fax: 86-754-8825 1499

E-mail: siui@siui.com

HONG KONG OFFICE:

Shantou Institute of Ultrasonic Instruments (HK) Co., Ltd.

Add: Room 2101, Tung Chiu Commercial Center

193 Lockhart Road, Wanchai, Hong Kong

Tel: 852-2891 6722 Fax: 852-2891 6723

www.siui.com

All rights reserved to SIUI 2016



SIUI



SIUI



SIUI

Apogee 3500 Elite/EN160331V1.0